

LED-4700 Battery Operated Headlight for Surgery

Soft Top
.....



Introducing our new product, the LED 4700! Designed with a keen focus on exceptional lighting performance, our latest headlight offers an unrivaled level of brightness and clarity while maintaining a lightweight and comfortable design. We are pleased to present this economic model with a variety of options that grants the ability to fine-tune the headlight to your specific needs.

Outshining the Competition in Brightness, Clarity and Comfort

- Advanced LED technology for exceptional brightness
- Light Intensity: 70,000 Lux (measured at a 16 in/40 cm standard working distance)
- Lightweight and ergonomic design for maximum comfort
- Premium Batteries built with a run time of 3 hours, 15 minutes at full intensity, 6.5 hours at half intensity
- Specialized Neoprene Wrap that is absorbent, breathable and easy to clean

Hard Top
.....




Our commitment to excellence shines through with the NEW LED-4700 Battery Operated Headlight. With this remarkable product, we have successfully combined superior quality with affordability, ensuring that a high value surgical headlight is accessible to all. Designed with the user in mind, the LED-4700 delivers a brilliant and bright light spot with excellent color representation. It is the perfect embodiment of our dedication to creating lighting solutions that empower surgeons while remaining comfortable and convenient to use.

Additional Specs:

- Color Temperature: 4,500° K for accurate color representation of tissue
- Typical LED life: > 50,000 hours
- Variable intensity control: 0 to 70,000 LUX (measured at a 16 in/40 cm standard working distance)
- Battery Recharge time 3.5 hours
- 3 year warranty for the LED headlight, 1 year for the batteries
- Audible warning when battery charge is low

System includes:

- 10 watt LED Headlight with choice of a Hard Headband or Soft Headband
- Battery Pack with two Standard Rechargeable Batteries
- Single Bay or a Dual Bay Battery Charger Option
- Reusable, autoclavable Joysticks for positioning light spot
- Carrying Bag for storage

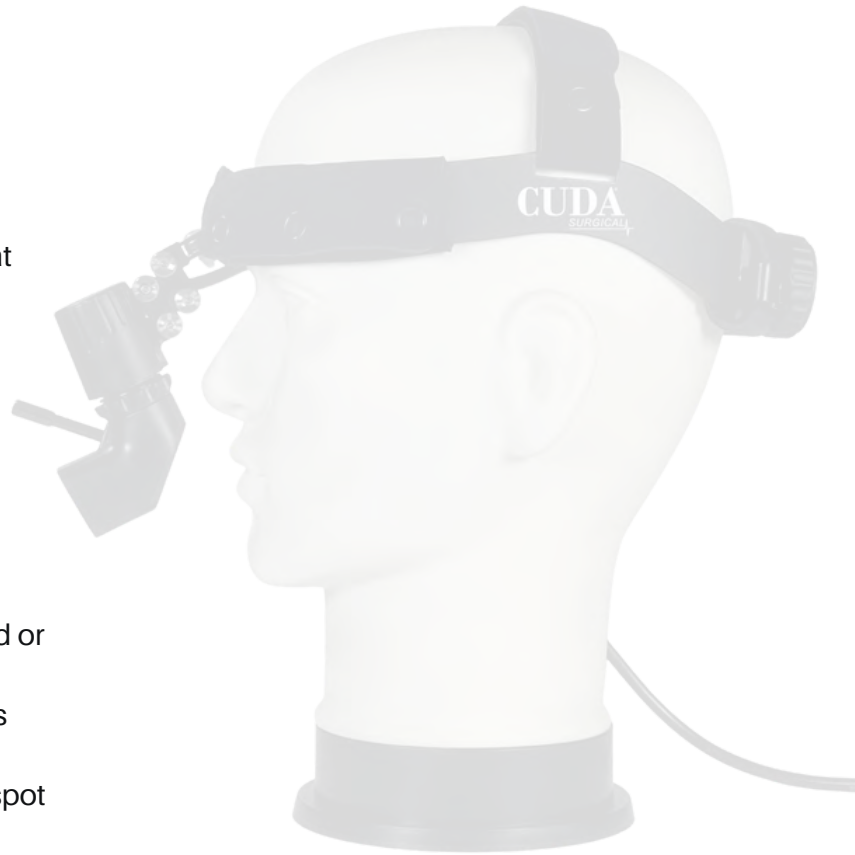

Options available:

- Choice of Hard or Soft Top Neoprene Wrapped Headband
- Single or Dual Bay Battery Charger
- 6 hour Extended Battery Option

Visit Our Website to Learn More:



www.CUDASurgical.com

Designed, Developed and Manufactured at



6018 Bowdendale Ave. Jacksonville, FL 32216, USA
Phone 877-814-2237 Fax: 904-733-4832

LIT 400

